



Vitalograph asma-1 range

the difference is clear

Vitalograph asma-1 range

Electronic asthma monitors

Vitalograph[®]
Your respiratory partner

Vitalograph asma-1

Electronic Asthma Monitor

- ❖ Simple to use
- ❖ Electronic record – no need for record cards
- ❖ Automatically assesses test quality
- ❖ More accurate than a mechanical meter
- ❖ Patient can review past results at home
- ❖ Measures PEF and FEV₁
- ❖ FEV₁ has less diurnal variability and is less effort dependent
- ❖ Can display % of personal best for both PEF and FEV₁
- ❖ PEF and FEV₁ zones can be personalised
- ❖ Automatically stores best values
- ❖ Large, easy to read display
- ❖ Stores 600 test sessions
- ❖ Suitable for use in clinic with SafeTway® mouthpieces
- ❖ Uses long life AAA batteries
- ❖ Easy to clean mouthpiece and flowhead
- ❖ Language independent
- ❖ ISO 23747:2007 Peak Flow Meters
ATS/ERS: 2005 Spirometry Guidelines



The asma-1® comes complete with:

- ❖ Plastic mouthpiece
- ❖ Carry pouch
- ❖ Guide for the healthcare professional
- ❖ 2 AAA batteries

The asma-1 is very easy to use, simply:

- ❖ Blow for > 1 second
- ❖ Repeat blow 3 times
- ❖ Press 'Enter' to view best result
- ❖ Transmit result (on comms versions)

Vitalograph asma-1 usb

For the home monitoring of asthma with electronic and hard copy reports

The asma-1 usb links to your PC to produce PDF reports for archiving or printing. These reports include the following:

- ❖ Subject ID and name
- ❖ Age, height, gender and optional weight field
- ❖ Best PEF and FEV₁ fields automatically populated
- ❖ BMI calculated (from height and weight)
- ❖ Personal best FEV₁, PEF and zones as set in asma-1
- ❖ Indicates blow quality, number of blows and diurnal variation
- ❖ Device ID automatically transmitted with the data
- ❖ Ability to add comments to the report
- ❖ Trend graphs of FEV₁ and PEF
- ❖ Optionally separate graphs of am and pm sessions
- ❖ Information presented in tabular and graphical form
- ❖ Can automatically name the PDF file
- ❖ Time/date stamp for each session



Additional Report Information

The following details are not stored on the model 4000 device. If you would like any of these details to be displayed on the generated report, please enter them now. Otherwise, they shall not appear on the report.


Demographics		Clinic	
First Name: John	Initial: EE	Date of Birth: dd/mm/yyyy: 16/09/1975	Clinic ID: 123456789012345
Last Name: Doe	Weight (kg): 75	Height (cm): 175	Clinic: _____
Subject ID: 123456789012345	Gender: Male	Unit Measure: <input type="radio"/> Imperial / US <input checked="" type="radio"/> Metric	Report Type: <input checked="" type="radio"/> All <input type="radio"/> Summary
Comments: _____		Graph AM/PM: <input type="radio"/> Together <input checked="" type="radio"/> Separate	

Printer Friendly Save Information



The asma-1 usb comes complete with:

- ❖ Reports Software & API/document CDs
- ❖ 4 disposable SafeTway mouthpieces and noseclips
- ❖ Plastic mouthpiece and carry pouch
- ❖ Guide for the healthcare professional
- ❖ USB cable and 2 AAA batteries

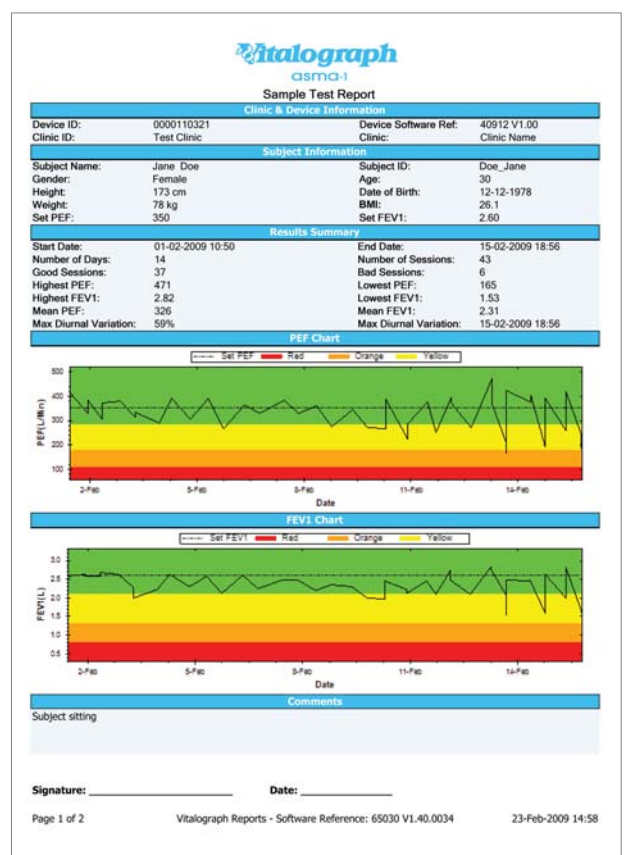

Technical Specification:

Product: Vitalograph asma-1
Vitalograph asma-1 usb

Model Number: 4000

Parameters Displayed: PEF & FEV₁

Memory: 600 PEF and FEV₁

Quality of Blow Indicator: Yes

Accuracy: FEV₁ better than +/-3%
PEF better than +/-5%

Range Volume: FEV₁ 0 - 9.99 L BTPS

Range Flow: PEF 25 – 840 L/min

Sensor: Stator-rotor mass flow meter

Flow Impedance: Better than 0.15 kPa/L/s
at 14 L/s

Power Supply: AAA batteries

Display: Custom liquid crystal display

Size: 113 x 63 x 48mm

Weight: 55g

Performance Standards:

ISO 23747:2007

ATS/ERS Guidelines: 2005

Operating Temperature: 17 - 37 °C

Management Zones: Factory set 3 zones
Green/Yellow threshold 80% Yellow/Red
threshold 50% Zones may be changed by the
user and number of zones increased to 4

Safety Standards: IEC 60601-1:2005

Medical Safety Standard: Medical Devices
Directive 93/42/EEC

Designed & Manufactured Under:

ISO 13485:2003, FDA 21CFR820, CMDR & JPAL

Ordering Info:

Catalogue No.

40000	Vitalograph asma-1
40400	Vitalograph asma-1 usb
40050	Vitalograph asma-1 s
40300	Vitalograph asma-1 bt
40200	Vitalograph copd-6
40550	Vitalograph lung monitor
40500	Vitalograph lung age
20242	SafeTway® mouthpieces (200)
20980	Mini SafeTway® mouthpieces (50)
28350	BVF (50)
20303	Disposable noseclips (200)
36020	3-L Precision syringe

4000 SERIES

	PEF	FEV ₁	FEV ₆	Ratio	Lung Age	Memory	Personal Best	GOLD	PDF Reports
asma-1	●	●	○	○	○	600	●	○	○
asma-1 (usb, s, bt)	●	●	○	○	○	600	●	○	●
copd-6™	○	●	●	●	●	○	○	●	○
copd-6 (usb, s, bt)	○	●	●	●	●	○	○	●	●
lung monitor	○	●	●	●	○	200	●	○	○
lung monitor (usb, s, bt)	○	●	●	●	○	200	●	○	●
lung age	○	●	○	○	●	○	○	○	○

● Available ○ Not available

usb - USB bt - Bluetooth s - Serial

Vitalograph®

Your respiratory partner

UK Sales
Vitalograph Ltd.
Maids Moreton Buckingham
MK18 1SW England
Tel: +44 1280 827110
Fax: +44 1280 823302
e-mail: sales@vitalograph.co.uk
www.vitalograph.co.uk

Export Sales
Vitalograph Ltd.
Maids Moreton Buckingham
MK18 1SW England
Tel: +44 1280 827120
Fax: +44 1280 823302
e-mail: sales@vitalograph.co.uk
www.vitalograph.eu

Vitalograph GmbH
Jacobsenweg 12
22525 Hamburg
Germany
Tel: (040) 547391-0
Fax: (040) 547391-40
e-mail: info@vitalograph.de
www.vitalograph.de

Vitalograph Inc.
13310 West 99th Street
Lenexa
Kansas 66215 USA
Tel: (913) 888 4221
Fax: (913) 888 4259
e-mail: vitcs@vitalograph.com
www.vitalograph.com

Vitalograph (Ireland) Ltd.
Gort Road Business Park
Ennis
Co Clare Ireland
Tel: (065) 6864100
Fax: (065) 6828289
e-mail: sales@vitalograph.ie
www.vitalograph.ie